



If you are now using an active threaded post...

Consider the advantages of ParaPostXT:
XTrA retention with XTrA safety!

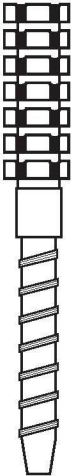

1. Rounded edges of the head minimize stress points that could lead to core failure.
2. Undercut, slotted head securely locks-on glass-ionomer and composite core materials firmly and easily.
3. Broad, flat shoulder stop provides security against over-insertion and the build-up of apical stress.
4. Patented, low-profile threads cut a clear path through the dentin, minimizing insertional stress. The interruptions permit venting of excess cement and dentin debris to relieve hydrostatic pressure.
5. Threads are located only in the stronger coronal section where root canal walls are thick. Apical trimming of post will not affect mechanical retention.
6. Passive, raised diamond-shape retention pattern in the apical section where root canal walls are more sensitive to fracture, provides additional retention and greater resistance against rotation.
7. Parallel sides distribute functional stress equally along the entire length of the post.

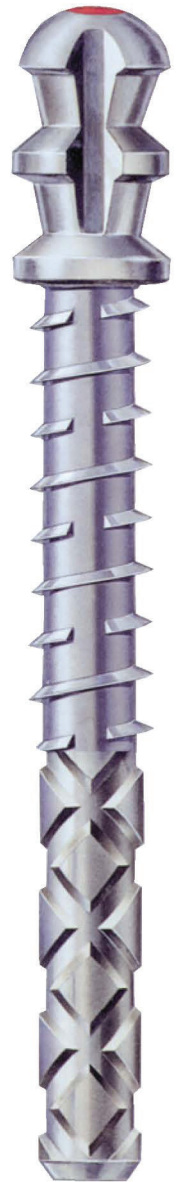

9062300 - Intro Kit

Refill

10/Pkg.

9062302 - Size 3 Brown

9062304 - Size 4 Yellow

9062306 - Size 4.5 Blue

9062308 - Size 5 Red

9062310 - Size 5.5 Purple

9062312 - Size 6 Black

Bulk Refill

30/Pkg.

9062313 - Size 3 Brown

9062314 - Size 4 Yellow

9062315 - Size 4.5 Blue

9062316 - Size 5 Red



DARBY DENTAL SUPPLY, LLC

To Reorder Call 800.645.2310 or Visit Darby.com
Contact your Account Manager for more information on ordering this product



XT

XT