


The All-New  
**A+Wedge**<sup>™</sup>  
Controls Gingival Bleeding




Insert

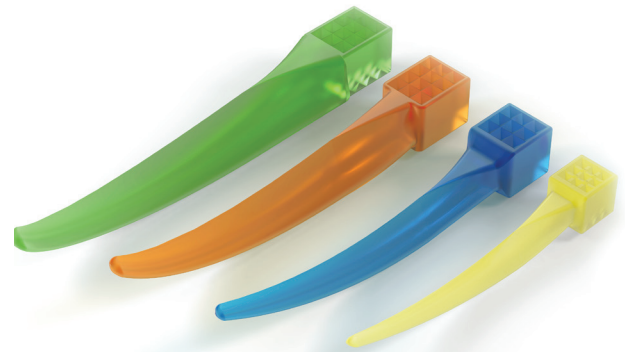

Wait & Restore



Remove & Rinse

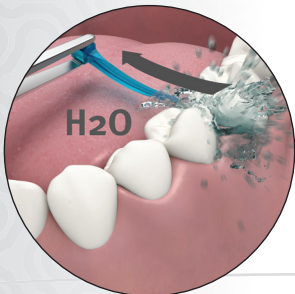

**The all-new A+ Wedge controls gingival bleeding** making restorations cleaner and quicker. In addition to adapting naturally to the contours of the teeth and mimicking interproximal anatomy, the A+Wedges are coated with aluminum sulfate that has an astringent effect on the gingiva. The wedges are designed to prevent back-out and are long enough to allow for proper matrix ring placement.

- Controls gingival bleeding
- Improves clarity of restoration area
- Eliminates need to apply an astringent separately

"The Garrison A+Wedge minimized or eliminated bleeding. It's particularly beneficial on patients with slight to mild gingivitis."

**Joe Arbutante, DDS**  
Garrison Advisory Board Member



**Garrison**  
Dental Solutions

# A+Wedge™

Controls Gingival Bleeding

More than your typical wedge, it's a wedge plus an astringent


**Stays put.** Back-out is prevented by the textured wedge body and gradual side taper.

**Longer is better.**


Matrix ring placement is easier, wedge back out is prevented, and fewer wedge sizes are needed with A+Wedge's more gradual increase in width.

**Snug fit.**

Adapts naturally to the contours of the teeth with the pliant plastic material.

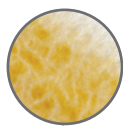

**Keeping it dry.** A+Wedge's astringent coating prevents bleeding while the continuous curve design reduces trauma. It's a must for all bonding procedures.

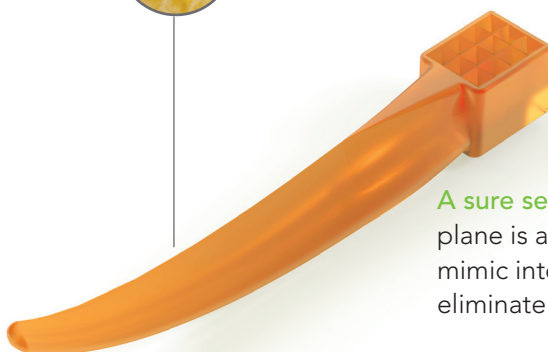

**A sure seal.** Every angle and plane is a continuous curve to mimic interproximal anatomy and eliminate flat spots.



**Gentle touch.** The ski-like tip prevents the wedge from digging into the papilla during placement, and the soft radius edges prevent trauma.



## Ordering Guide

### A+Wedge™

Kit	Description	Qty.
8390241 (GWAK4)	A+Wedge Kit	50/ea. size
<b>Refills</b>		
8390243 (GWAYL)	X-Small Yellow Refill	100
8390245 (GWABL)	Small Blue Refill	100
8390247 (GWAOR)	Medium Orange Refill	100
8390249 (GWAGR)	Large Green Refill	100

To Reorder Call 800.645.2310 or Visit [darby.com](http://darby.com)