

AFFINIS®

PERFECT IMPRESSIONS

BLACK EDITION

2 x 75 ml = 4 IMPRESSIONS

AFFINIS® heavy body Single Package

2 x 75 ml heavy body BLACK EDITION **9068212**
8 x Mixing Tips **(6001 9554)**

including
PLUNGER

AFFINIS® Trial Kit

2 x 75 ml heavy body BLACK EDITION
2 x 50 ml PRECIOUS light body/regular body **9068214**
1 x Plunger 75 ml **(6001 9555)**

AFFINIS® System 360 heavy body Starter Kit

1 x 380 ml heavy body BLACK EDITION **9068216**
1 x Fixation Ring
10 x Dynamic Mixing Tips - Yellow **(6001 9774)**

AFFINIS® System 360 heavy body Refill

2 x 380 ml heavy body BLACK EDITION **9068218**
1 x Fixation Ring **(6001 9776)**

AFFINIS® PRECIOUS wash material Gold

2 x 50 ml regular body GOLD **9068185 (6776)**
12 x Mixing Tips
4 x 25 ml regular body GOLD **9068187 (6780)**
20 Mixing Tips | 20 Oral Tips

AFFINIS® PRECIOUS wash material Silver

2 x 50 ml light body SILVER **9068184 (6775)**
12 x Mixing Tips
4 x 25 ml light body SILVER **9068186 (6779)**
20 Mixing Tips | 20 Oral Tips



© 2015 Coltène/Whaledent AG



DARBY DENTAL SUPPLY, LLC

To Reorder Call 800.645.2310
or Visit darby.com

COLTENE
coltene.com

P/N 31271B

AFFINIS[®]
PERFECT IMPRESSIONS



BLACK
WITH LOVE FOR DETAILS

 **COLTENE**



OUTSTANDING READABILITY


The new black tray material AFFINIS heavy body BLACK EDITION provides excellent readability of the impression when used with gold or silver colored AFFINIS PRECIOUS.

The unique color combination displays details even more precisely and supports the qualitative assessment of the impression result.

Visible precision with love for details.

The unique silver and gold pigmentation drastically reduces the light scattering effect and the readability is intensified by maximum color contrast.


AFFINIS® PRECIOUS

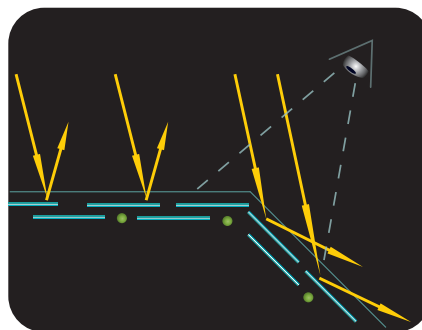

Filler ●
AFFINIS® PRECIOUS pigments —


The lamellar structure of the mirror like filler reduces light scattering of the surface which results in a better contrast and readability.

Standard impression material

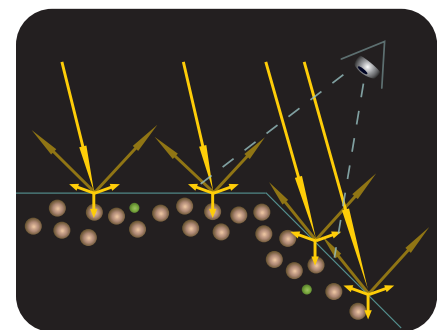

Filler ●
Standard pigments ●

The irregular structure of the standard fillers leads to undirected light scattering at the surface which results in less contrast and inferior readability.

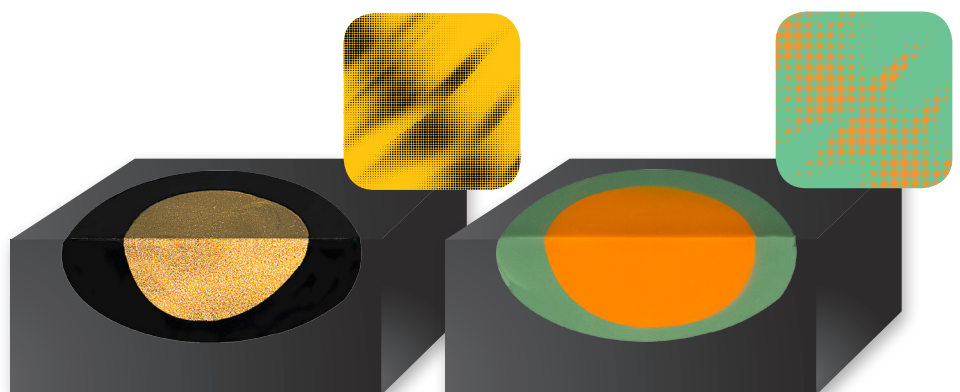

SUPERIOR FLOW PROPERTIES


AFFINIS® PRECIOUS wash materials have a contact angle of only 10°. Their superior flow properties allow the impression material to easily capture all critical surface

detail of the preparations in a wet environment. The unique surface affinity prevents air entrapments and margin distortion.

For perfect impressions.




Initial surface affinity on moist dentin

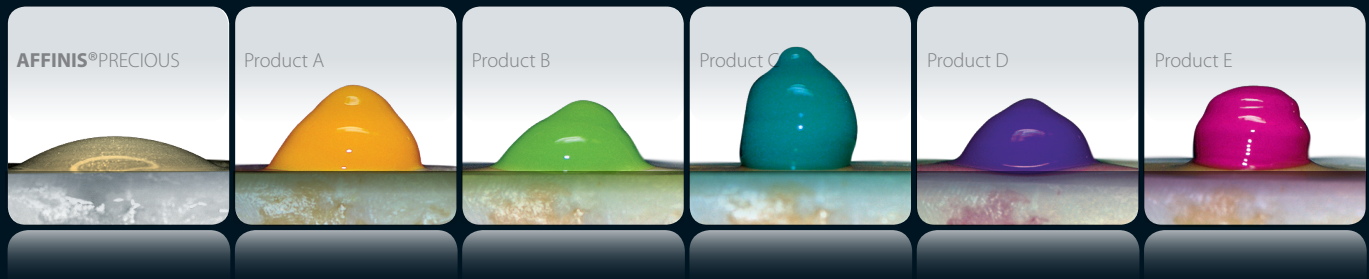

Defined amount of material on moist dentin surface after 10 s.

AFFINIS® PRECIOUS was developed to produce accurate results in the mouth and in the lab. The wettability allows the wet gypsum to flow more easily into the impression, producing more accurate stone models. An accurate model is the key for the lab technician to create properly fitting restorations.

AFFINIS® PRECIOUS regular body

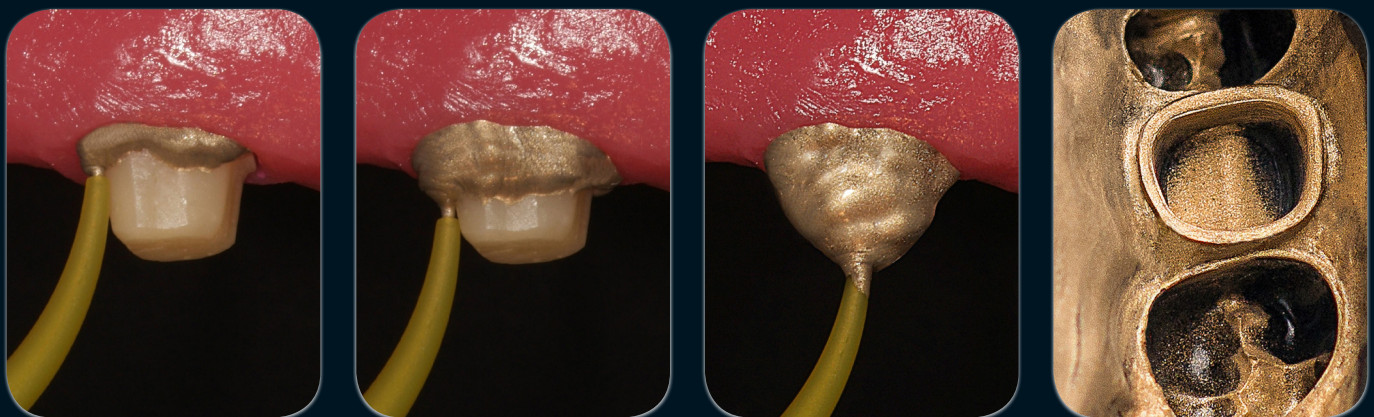

AFFINIS® PRECIOUS flows directly into the sulcus preventing air entrapments.

Thixotropic properties prevent slumping.

The self-contouring consistency and the highly active surfactant agent produces an innovative affinity to the moist tooth structure.

AFFINIS® PRECIOUS wash materials produce accurate, void free impressions.